



4.05 Cabinet Oxidising Agent Storage 100L

Self closing mechanism Easily replaced should door be damaged

Full instructions on usage of the cabinet attached to the inside of the door

Round corners on door to eliminate cuts

Welded shelf brackets that fully interlock with shelves. 40mm gap between brackets

Danger alert sign indicating specific hazardous substances stored in cabinet

Fully adjustable galvanised shelves perforated for any leakages and to permit free movement of air. Easily replaced if damaged

Ventilation port with built in flame arrestor

Self closing non latching door held shut by a magnetic lock which releases in the event of an internal buildup of pressure

Easy grip handle for opening of cabinet. Lock keyed to 92268


Ventilation port with built in flame arrestor

Continuous piano hinge ensuring smooth closure of door

At 800mm high designed to be installed under a standard laboratory bench top

Fully welded liquid tight sump capable of containing the entire contents of the cabinet

Baked enamel lead free paint inside and out eliminating corrosion issues.

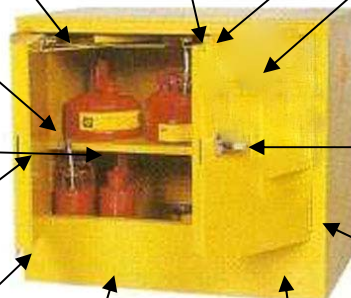

Model	Capacity (L)	Shelves	External Height (mm)	External Width (mm)	External Depth (mm)	Internal Height (mm)	Internal Width (mm)	Internal Depth (mm)	Weight (kg)
BCOAS 100L	100	1	800	920	615	520	835	530	100

Australian Standard

All cabinets comply with AS4326- 1995 as follows.

1. Each cabinet is marked with the name and address of the manufacturer
2. The maximum storage capacity.
3. A sign bearing the words "NO SMOKING, NO IGNITION SOURCES WITHIN 3M" in lettering at least 50mm high
4. Maximum quantity for storage of PG1 substances is 50kg/lt

Options

Full range of spare parts including additional shelves, closing mechanisms etc. available.

Storage

Suitable for safe storage of oxidizing agents as defined in the ADG code. Cabinets used for storage of oxidizing agents shall not be used for storage of other materials. Where more than one oxidizing agent is stored in a cabinet, the oxidizing agents shall be compatible.

Made in Australia

Manufactured by Brodie in Australia using Australian products to suit Australian conditions.



Brodie

Head Office 42 Lisbon Street Villawood 2163
 Phone 61 2 9722 4600
 Fax 61 2 9722 4666
 Email sales@aafab.com.au
 Website www.brodie.net.au